Accessory Kit Instructions for Use

- EN English
DE Deutsch
FR Français
ES Español
NO Norsk
SV Svenska
AR عربي

Aerogen Ltd.
Galway Business Park,
Dangan, Galway,
Ireland, H91 HE94
Tel: +353 91 540 400
E-mail: info@aerogen.com
www.aerogen.com

HCmed Innovations Co., Ltd.
Rm.B, 10F, No.319, Sec.2,
Dunhua S. Rd.,
Taipei City, 10669, Taiwan
Tel: +886-2-2732-6596
E-mail: info@hcmmed-inno.com
HCM-86C

EC REP Medical Device Safety Service GmbH,
Schiffgraben 41, 30175 Hannover,
Germany
CH REP Medical Device Safety Service CH GmbH,
Laurenzenvorstadt 61, 5000 Aarau,
Switzerland

UK Responsible Person:
MDSS UKRP Ltd., 6 Wilmslow Road,
Rusholme, Manchester M14 5TP,
United Kingdom

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
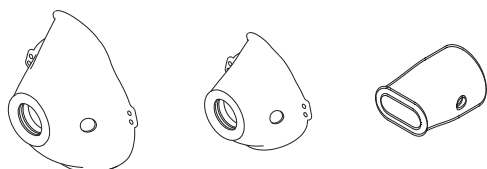
Lagerbedingungen: -25~70 °C, max. 93 % relative Luftfeuchtigkeit, 700~1060 hPa

Accessory Kit Instructions for Use

EN

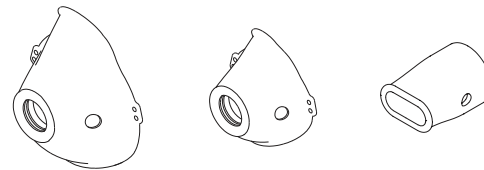Product Name	Aerogen® Go Mask (Adult)	Aerogen® Go Mask (Paediatric)	Aerogen® Go Mouthpiece
Product No.	AG-AG8021-02	AG-AG8022-02	AG-AG8031-02
Model No.	DP-M001	DP-M002	DP-P001

Caution:

This kit is designed for a single patient. Do not share with other patients. This instruction is for general users. Babies, children and anyone who require assistance or supervision must be observed constantly by an adult or medical professional when using the accessory for inhalation therapy.

Indication: To assist inhalation therapy with Aerogen® Go Vibrating Mesh Nebulizer.

Contraindication: There is no contraindication for this kit. However, failure to comply with the Instructions for Use can lead to weak or erratic aerosol stream, and can impair effectiveness.

Instruction:

- Connect either the mouthpiece or the mask to the Aerogen® Go medication reservoir.
- Make sure the accessory is well connected to the Aerogen® Go medication reservoir before use.

Cleaning and Disinfection: For more details, please refer to pages 10 of Aerogen® Go Vibrating Mesh Nebulizer Instruction Manual.

Expected lifetime: The Accessory Kit is designed to provide one year use with 30 minutes operation per day. The life cycle of the components may be shortened if the Aerogen® Go Vibrating Mesh Nebulizer is used for longer durations or extended if the Aerogen® Go Vibrating Mesh Nebulizer is used for shorter durations.

Operating condition: 5°C~40°C, 15~93%R.H. (non-condensing), 700-1060hPa

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30-2175 Rev 01 © 2025 Aerogen Ltd.